REVANESSE® CONTOUR™ XL

Instructions for Use
Gebrauchsanweisung
Mode d'emploi
Istruzioni per l'uso
Instrucciones de uso
Instruções de Uso
Инструкция по применению

SYMBOL SYMBOL SIMBOLE SIMBOLO
SIMBOLO SIMBOLOGIA SIMBOLIT:

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Date of Manufacture
Date de fabrication

Keep dry
Garder au sec

Expiration date
Date d'expiration


Manufacturer
Fabricant

ASSEMBLY OF NEEDLE TO SYRINGE:

1. Open bag and verify the kit contents (1×10 mL syringe, 1×plunger rod, 1× finger grip, 1×instructions for use).
2. Select an appropriate cannula for the injection of dermal filler (sold separately).
3. Hold the 10 mL syringe around the barrel and with the other hand, snap the finger grip into place with its logo facing up.
4. Continue to hold syringe and align the plunger rod with the screw facing down.
5. Place the screw with gentle pressure on to the rubber stopper and affix it with a clockwise motion (stop when you feel resistance).
6. Unscrew the tip cap of the syringe and affix the cannula by holding the syringe luer-lock and affix it in a clockwise motion.

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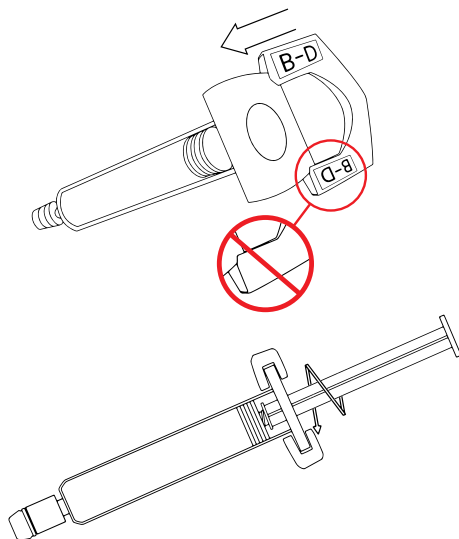COMPOSITION

Cross-linked hyaluronic acid (High Viscosity) 25mg/ml
 In phosphate buffered saline
 (Cross-linked with Butandiol-diglycidylether (BDDE))

DESCRIPTION

Revanesse® Contour™ XL is a colorless, odorless, transparent and aqueous gel of synthetic origin. The gel is stored in a pre-filled disposable syringe. Each box contains one 10ml syringe of Revanesse® Contour™ XL.

APPLICATION RANGE / INDICATIONS

Revanesse® Contour™ XL is a cross linked hyaluronic acid gel that is indicated for the treatment of medium to deep facial rhytides and lipatrophy by injection into the middle part of the dermis layer. Implantation life is dependent on depth and location of injection, and averages 6-12 months.

ANTICIPATED SIDE EFFECTS

Physicians must inform patients that with every injection of Revanesse® Contour™ XL there are potential adverse reactions that may be delayed or occur immediately after the injection. These include but are not limited to:

- Injection-related reactions might occur, such as transient erythema, swelling, pain, itching, discoloration or tenderness at the injection site. These reactions may last for one week.
- Nodules or induration are also possible at the injection site.
- Poor product performance due to improper injection technique.
- Glabellar necrosis, abscess formation, granulomas and hypersensitivity have all been reported with injections of hyaluronic acid products. It is important for physicians to take these reactions into account on a case by case basis. Reactions thought to be of hypersensitivity in nature have been reported in less than one in every 1500 treatments. These have consisted of prolonged erythema, swelling and induration at the implant site.

These reactions have started either shortly after injection or after a delay of 2 – 4 weeks and have been described as mild or moderate, with an average duration of 2 weeks. Typically, this reaction is self-limiting and resolves spontaneously with time. However, it is imperative that patients with hypersensitivity type reactions contact their physician immediately for assessment. Patients with multiple allergic reactions should be excluded from the treatment.

CONTRAINDICATIONS

- Do not inject Revanesse® Contour™ XL into eye contours (into the eye circle or eyelids).
 - Pregnant women or lactating women should not be treated with Revanesse® Contour™ XL.
 - Revanesse® Contour™ XL is only intended for intradermal use and must not be injected into blood vessels. This may occlude and could cause an embolism.
 - Patients who develop hypertrophic scarring should not be treated with Revanesse® Contour™ XL.
 - Patients that have a known history of streptococcal disease.
 - Contains trace amounts of gram positive bacterial proteins, and is contraindicated for patients with a history of allergic reactions to this material.
 - Never use Revanesse® Contour™ XL in conjunction with a laser, intense pulsed light, chemical peeling or dermabrasion treatments.
 - People under the age of 18 should not be treated with Revanesse® Contour™ XL.
 - Patients with acne and / or other inflammatory diseases of the skin should not be treated with Revanesse® Contour™ XL.
 - Patients with unattainable expectations.
 - Patients with auto-immune disorders or under immunotherapy.
 - Patients with multiple severe allergies.
 - Patients with acute or chronic skin disease in or near the injection sites.
 - Coagulation defects or under anti-coagulation therapy.
 - Patients with sensitivity to hyaluronic acid.
- It is imperative that patients with adverse inflammatory reactions that persist for more than one week report this immediately to their physician. These conditions should be treated as appropriate (ie: corticosteroids or antibiotics). All other types of adverse reactions should be reported directly to the authorized distributor of the Revanesse® family of products and / or to ProLium Medical Technologies Inc. directly.

ADMINISTRATION & DOSAGE

- Revanesse® Contour™ XL should only be injected by or under the direct supervision of qualified physicians who have been trained on the proper injection technique for filling facial wrinkles.
- Before patients are treated they should be informed of the indications of the device as well as its contraindications and potential undesirable side effects.
- The area to be treated must be thoroughly disinfected. Be sure to inject only under sterile conditions.
- Inject the product slowly and apply the least amount of pressure necessary.
- Revanesse® Contour™ XL and needles packaged with it are for single use only. Do not re-use. If re-used, there is a risk of infection or transmission of blood born diseases.
- Keep the product at room temperature for 30 minutes prior to injection.
- If the skin turns a white color (blanching), the injection should be stopped immediately and the area should be massaged until the skin returns to its normal color.
- Before injecting, press on the plunger of the syringe until a small drop is visible at the tip of the needle.

PRECAUTIONS

- Revanesse® Contour™ XL should not be injected into an area that already contains another filler product as there is no available clinical data on possible reactions.
- Revanesse® Contour™ XL should not be injected into an area where there is a permanent filler or implant.
- Hyaluronic acid products have a known incompatibility with quaternary ammonium salts such as benzalkonium chloride. Please ensure that Revanesse® Contour™ XL never comes into contact with this substance or medical instrumentation that has come into contact with this substance.
- Revanesse® Contour™ XL should never be used for breast enlargement, or for implantation into bone, tendon, ligament or muscle.
- Avoid touching the treated area for 12 hours after injection and avoid prolonged exposure to sunlight, UV, as well as extreme cold and heat.
- Until the initial swelling and redness have resolved, do not expose the treated area to intense heat (e.g. solarium and sunbathing) or extreme cold.
- If you have previously suffered from facial cold sores, there is a risk that the needle punctures could contribute to another eruption of cold sores.
- If you are taking aspirin, non-steroidal anti-inflammatory medications, St. John's Wort, High doses of Vitamin E supplements or any similar medications prior to treatment, be aware that these may increase bruising and bleeding at the injection site.
- The safety for use in patients under 18 years or over 65 years has not been established.
- Patients who are visibly ill, with bacterial or viral infections, influenza, or active fever should not be treated until a resolution of their symptoms.
- In order to minimize the risks of potential complications, this product should only be used by health care practitioners who have appropriate training, experience, and who are knowledgeable about the anatomy at and around the site of injection.
- Health care practitioners are encouraged to discuss all potential risks of soft tissue injection with their patients prior to treatment and ensure that patients are aware of signs and symptoms of potential complications.

WARNINGS

- Confirm that the seal on the box has not been broken and sterility has not been compromised. Confirm that the product has not expired. Product is for single use only; do not re-use. If re-used, there is a risk of infection or transmission of blood born diseases.
- Introduction of product into the vasculature may lead to embolization, occlusion of the vessels, ischemia, or infarction.
 - Rare but serious adverse events associated with the intravascular injection of soft tissue fillers in the face have been reported and include temporary or permanent vision impairment, blindness, cerebral ischemia or cerebral hemorrhage, leading to stroke, skin necrosis, and damage to underlying facial structures.
 - Immediately stop the injection if a patient exhibits any of the following symptoms, including changes in vision, signs of a stroke, blanching of the skin, or unusual pain during or shortly after the procedure.
 - Patients should receive prompt medical attention and possibly evaluation by an appropriate health care practitioner specialist should an intravascular injection occur.
- Physicians should take extra precautions when treating patients with auto immune disorder, HIV, AIDS or who are under immuno-therapy as these patients are at higher risk of adverse events.
- The Revanesse® family of products should not be used in areas that have high vascularity. Use in these areas such as the glabella and nose region has resulted in cases of vascular embolization and symptoms consistent with ocular vessel occlusion (i.e.: blindness).
- NOTE:** The correct injection technique is crucial to treatment success & patient satisfaction. Revanesse® Contour™ XL should only be injected by a practitioner qualified according to local laws and standards.
- The graduation on the syringe is not precise and should be used as a guide only. The amount of material to be injected is best determined by visual and tactile assessment by the user.

SHELF LIFE & STORAGE

Expiry is indicated on each individual package. Store between 2°-25° C, and protect from direct sun light and freezing.

MANUFACTURER

ProLium Medical Technologies, Inc.
 138 Industrial Parkway North, Aurora, ON
 L4G 4C3, Canada

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Prolenium Medical Technologies, Inc.